
# PRIXTON

MGF100

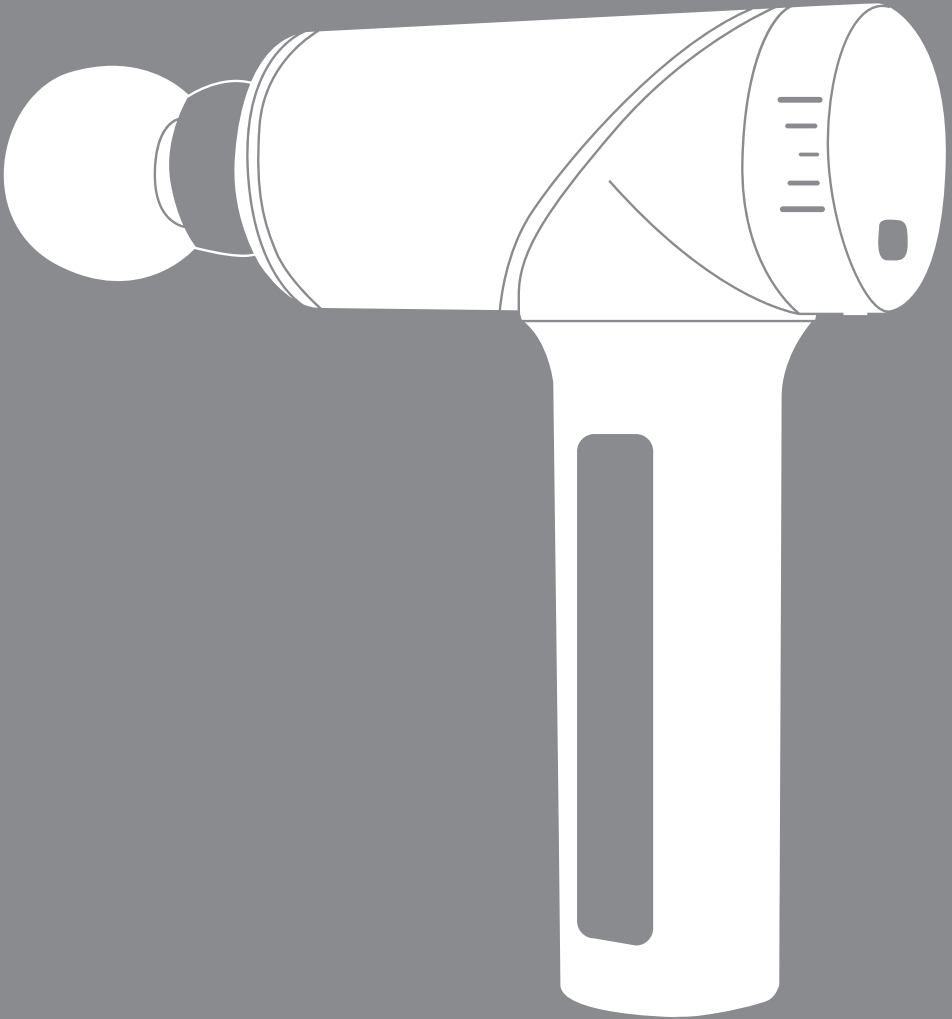**PISTOLA MASAJE**

**MASSAGE GUN**

---

## Declaración

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---

**Statement**

To reduce the risk of faults, electric shocks, injuries, death or damages caused by fire, please read carefully all the instructions, advice and warnings in this manual before using the product.

**\*Warnings and precautions**

1. Minors must not use this product if they have injuries. Check with your doctor before using this device if you have a muscular injury.
2. Do not use this product on your skin unless it is clean and dry. You can use it on clothes. Press and move the device gently for approximately 60 seconds.
3. Do not use the device above the neck and on parts of the body with little muscle mass or on bone. It may only be used on soft tissue free from pain and injuries.
4. Do not use the product on bruised skin or where there is pain or irritation.
5. Do not immerse the product in water.
6. Do not put the product near a flame or source of fire.
7. Do not let the product fall onto the floor.
8. Children aged under 13 should only use it under adult supervision.
9. Do not take the product apart in any way.
10. Keep fingers, hair, and other body parts away from the axle or massage head, because they can become entangled.

**\*Warnings and reminders**

1. Keep the product away from liquids or sources of heat. Keep the product dry.
2. After charging, unplug the charger.
3. Only use the heads that give the best results.
4. Use the device in accordance with the instructions.
5. Please use the original charger to charge the battery.
6. Check the condition of the device before each use.

**\* Important information**

This product must not be used in the following situations without prior consent from a doctor:

People with complications related with pregnancy and urination, neuropathology or retinal injury , people wearing pacemakers or who have been operated on recently, epilepsy or migraine sufferers, people with slipped disks, spondylolisthesis, dissection of the vertebral artery or rigidity in the articulation of the vertebral arch, who have recently been operated on to replace articulations or metal needle or steel plate.

The presence of these conditions does not mean that you cannot use the product. We suggest you consult your doctor before using it. Current research focuses on training with vibrations for specific illnesses. It will probably cure the contra-indications listed above. Experience shows that the combination of training with vibrations and long-term therapy is a viable form of treatment in many cases. It should be applied following the instructions or with supervision by a doctor, expert or physiotherapist.

## Product description

1. Powerful motor without brushes, with a silent operating technology.
2. Metal heat diffuser.
3. The personalized gentle heads provide a pleasant massage and reduce injuries to the body.
4. Recharge the Lithium -Ion battery. The battery cannot be removed. Do not extract. Do not remove by force.
5. Nine power levels and two massage modes to use as you prefer.
6. Wireless, light and comfortable. The handle is made of silicon gel, making it smooth and comfortable to use.
7. Charging dock. Please charge the device using the charger provided.

1. Fitness mode (FITNESS): 9 levels of constant speed. High speed and torque, suitable for fitness and professional athletes. Ideal for use as a warm-up before and after training.
2. Physiotherapy mode (MEDI): 2 variable speed levels. Increases and decreases the force and speed to treat a local area or therapy application.
3. ON/OFF button: Press for a few seconds to turn the gun on or off. Press once to increase speed. To switch between FITNESS and MEDI modes, double press twice.
4. Speed indicator: Levels from 1-9 in FITNESS mode and 1-2 in MEDI mode.
5. Battery indicator: Displays the charge level in percentage. When the light flashes, it indicates low battery level. NOTE: To preserve the life of the battery, do not allow it to discharge completely.

**With 6 interchangeable massage heads.**

**Spherical head**  
Use to massage large muscle groups, such as arms, back, thighs and calves

**Air cushioned head**  
Use for massaging painful areas. Used to reduce vibrations.

**Conical head**  
Can be used to treat deep tissue such as the meridian points, the palms of the hands, the soles of the feet, etc.

**Y-shaped head**  
Suitable for massaging the neck, the spine, Achilles' tendon, etc.

**Golden fingers**  
Used to relax different groups of muscles, such as the meridian points, palms of the hands, etc. To massage the neck, spine, Achilles' tendon, etc.

**Flat head**  
Suitable for relaxing and shaping various muscle groups.

## Charging

1. Before using the device for the first time, charge the device for several hours to obtain a full charge.
2. Please charge the device using the charger provided.
3. The LED lights of the device will flash regularly to indicate the battery level and that the device is charging.
4. If the indicator light of 25% is lit, there is between 0 and 25% of the battery power left. If the indicator light of 50% is lit, there is between 26 and 50% of the battery power left. If the indicator light of 75% is lit, there is between 51 and 75% of the battery power left. If the indicator light of 100% is lit, there is between 76 and 100% of the battery power left.
5. The battery may be charged at any time to any level. The device cannot be used when it is charging.
6. When the battery charge is below 25% (the LED light for 25%), it must be charged.
7. The average time required for a full charge is between 3 and 4 hours, and there are differences in speed and strength.

## Installing the massage head

Turn the device off, align the head with the mouth of the massage pistol. Insert the part slowly (in the direction indicated by the red arrow in image 1) and attach it as shown in image 2.

Image 1. Direction for installation

Image 2. Installation completed

## Maintenance, cleaning, storage and transport

1. Clean the device with a towel or kitchen roll and dry it with a smooth cloth.
2. When storing the product or taking it when travelling, put it in the packaging provided.
3. It must be stored in a dry, well-ventilated space, free from corrosive gases. It should be no less than 250mm from the floor.
4. If it is not used for a long period, it must be wrapped up and put away after cleaning it. The product must be charged every 6 months.
5. The box can be carried by car, train, boat or plane.
6. In transport, fasten the product to prevent strong vibrations, yawing and falls. It must be raised, handled gently and protected from the rain.

**Important:**

1. To avoid injuries, do not insert fingers or other objects into the gap between the massage head and the body.
2. Do not overuse the massage pistol.
3. Do not apply the pistol on bones.
4. Do not use this device without the consent of your doctor if you are pregnant, a child or teenager, if you are suffering complications, injured, wear a pacemaker or suffer from osteoporosis.
5. If you feel pain or feel unwell, please stop using the device at once.

**Warranty policy**

This product is guaranteed for 2 years from the date of purchase. When you have a problem that you can not solve, you must access the web [www.prixton.com](http://www.prixton.com) and click on the contact option to send us your assistance form.

**Simplified eu declaration of conformity**

We, La Trastienda Digital located in Avenida Cervantes 49, módulo 25, Basauri 48970, Bizkaia hereby declare that the following goods, imported and selling at European market:

| <b>BRAND</b> | <b>MODEL NAME</b> | <b>DESCRIPTION</b> | <b>CATEGORY</b> |
|--------------|-------------------|--------------------|-----------------|
| Prixton      | MGF100            | Massage Gun Fit    | Fitness         |

Are in compliance with the EC Council Directive of 2014/53/EU (EM DIRECTIVE) and harmonized standards and with the board of the 16th of April 2014.

The complete text of the UE Declaration of conformity is available anytime at this website address:  
[www.prixton.com](http://www.prixton.com)

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